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
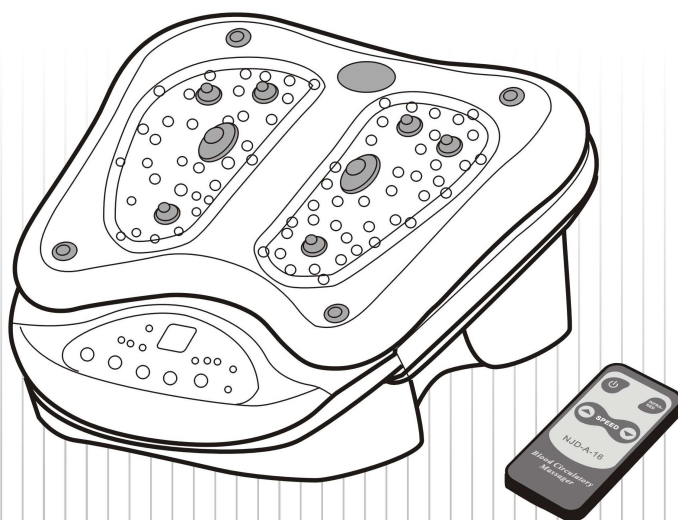
**DANDELION Massager Vibration Machine**
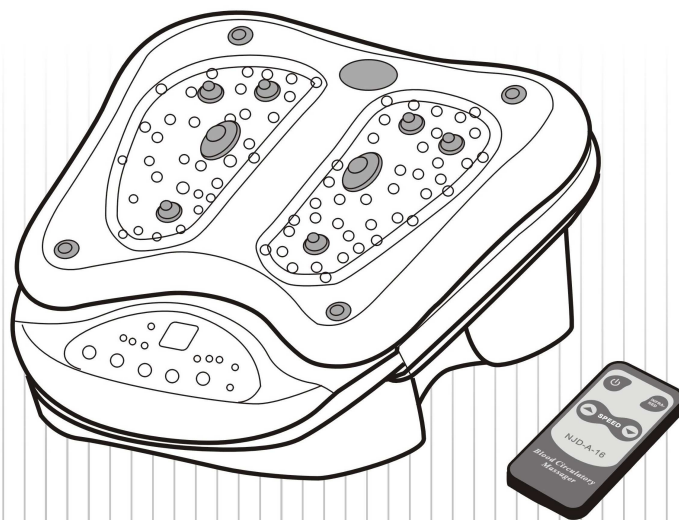
## **GENERAL SAFETY ATTENTION**

**DANGER:** to reduce the risk of electric shock:

- Always unplug the device immediately after using and before each cleaning. When plugged in, it should never be left unsupervised.
- Do not attempt to pick up a device that has fallen into water or any other liquid. Unplug it immediately.
- Do not use this device in humid or extremely dirty environments, as this may cause a power cut or result in incorrect operation of the machine.
- Do not use this device when taking a shower or bath.
- Do not place or store this device in any place from where it might fall or be pulled into a bath or sink. Refrain from placing it or dropping it in water or any other type of liquid.

**ATTENTION:** to reduce the risk of burns, fire, electric shock or bodily injury.

- Close supervision is necessary when this device is

used by, on or near children, disabled people or handicapped people.

- Use this device only for its intended purpose, as described in the operation manual. Use only the accessories supplied with this unit.

## **GENERAL SAFETY ATTENTION**

- Refrain from using any device where the cord or plug has been damaged, as well as any machine that is not working correctly, that has fallen or is damaged, or that has fallen into water. Send it to a qualified technician or a service centre for inspection and repair.
- Keep the cord away from hot surfaces.
- Never use the device while you are sleeping. Never fall asleep when you are using the masseur.
- Never drop the unit and never insert objects into its openings.
- Refrain from using the equipment when using aerosol products (sprays) or when oxygen is being administered.
- Refrain from using the equipment on or near the eyes or any other extremely sensitive area.

- Refrain from using the equipment under covers or pillows. This may produce excessive heat and cause a fire, electric shock or harm the body.
- Do not move this device by using the power supply cord and refrain from using the cord as a handle.
- Do not use the device outside.
- To unplug it, turn the power supply to the "OFF" position; then pull out the plug.
- It is forbidden for children to use this unit.
- Do not cover the device while it is in operation.
- Do not use this product for more than 20 minutes in any single session.
- This device is intended only for domestic use.

## **GENERAL SAFETY ATTENTION**

- If you have the slightest doubt about your health,  
please consult your doctor before using this product.  
The following users should obtain approval from

their doctor before using the machine:

- Patients having suffered a brain clot.
- Patients with a heart condition.
- Patients with a deformation of the vertebral column  
or receiving treatment for a fracture or fracture  
deformity.
- Pregnant women, menstruating women or people  
suffering from fever following an operation.
- People who are mentally impaired.
- Patients with acute illnesses or serious cancer.
- Infectious patients with skin diseases.
- Patients still recovering from an operation.
- People wearing a pacemaker should consult a doctor  
before use.
- Diabetics.
- If the supply cord is damaged, it must be replaced by  
the manufacturer, its service agent or similarly  
qualified persons in order to avoid a hazard.
- This appliance is not intended for use by persons

(including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.


## **GENERAL SAFETY ATTENTION**

- Prolonged use could lead to overheating of the product and shorten its useful lifetime. To avoid these inconveniences, stop the unit and let it cool down before reusing it.
- Never use this product directly on swollen or inflamed areas or areas with open wounds or skin rashes.
- This machine is not intended for professional use. Do not use this product in medical treatment areas.
- Do not use this product before going to sleep. The massage has a stimulating effect and may cause insomnia.
- This product should never be used by any person suffering from physical problems likely to hinder their capacity to use the controls.
- If the power supply cord should become damaged, please see a qualified technician for the necessary repairs or maintenance to avoid any risks.
- If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a

hazard.

- This product should be used under the feet of a sitting person.

## **PARTS NAME**

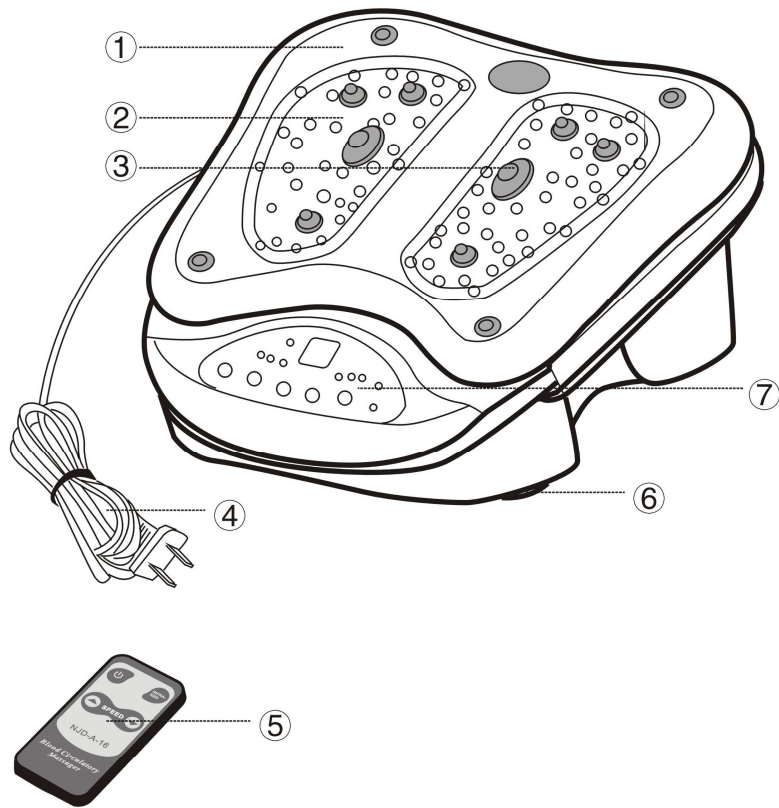①Vibrating plate

②Infrared light

③Massage points

④Cable

⑤Remote control

⑥Support Stand

⑦Control panel

## OPERATION MANUAL

1. Place the device on a flat surface and plug it into the electric socket.

2. Press the "ON/OFF" button; the device will begin working with the pilot light being illuminated. The LCD screen displays the value 6; this means that the intensity of the vibrations is in position 6. The initial setting is in manual mode: vibration speed 6, timer at 15 minutes and infrared heat treatment at speed 3.

2. Press ▲ to increase the rate of the vibration and ▼ to reduce this rate. Every time it is pressed, there will be either an increase or decrease in the level. At the same time, the screen will display the number corresponding to the change in rate. The control range for the speed goes from 0 to 12, with 0 meaning the massage will end.

3. Press "Time" to select the duration of the operation. Every time you press this you can move successively through the modes 10 minutes 15 minutes 20 minutes and 10 minutes. The initial setting is 15 minutes.

4. Press on "INFRARED" to change the thermotherapy. The operation range for the infrared light varies between H0 and H5, where H0 indicates the absence of heat. The initial setting is at H3. Each time this is pressed you can successively move through modes H4-H5-H0-H1-H2-H3-H4. Every time you press "INFRARED" you are able to read the heat level

on the LCD screen. The screen will again display the speed level after a period of 8 seconds each time it is pressed.

## **OPERATION MANUAL**

5. Press "AUTO" to select the automatic mode.

This unit has two programmes: circulation I and circulation II. Each time you press "AUTO" you can successively move through the modes Manual Circulation I Circulation II Manual.

(1) Circulation I: Possibility of adjusting the vibration speed between 2 and 12. Possibility of adjusting the infrared heat between 5 and 1. The time between each change in mode is 8 seconds.


(2) Circulation II: Possibility of successively adjusting the vibration speed between 2-8-4-6-12-7-9-5. Possibility of adjusting the infrared heat: 5-5-4-4-1-1-3-3. The time between each change in mode is 4 seconds.

6. Select the appropriate time, temperature and speed depending on your state of health.

We recommend using the device between 1 and 3 minutes for every initial use in a given place. Following this, change the time, temperature and speed according to your wishes. It is vital that you select the functions beginning with a low level of intensity.

7. Press "ON/OFF" to stop the power supply then pull the plug out..

## REMOTE CONTROL: TO BE USED IN MANUAL MODE
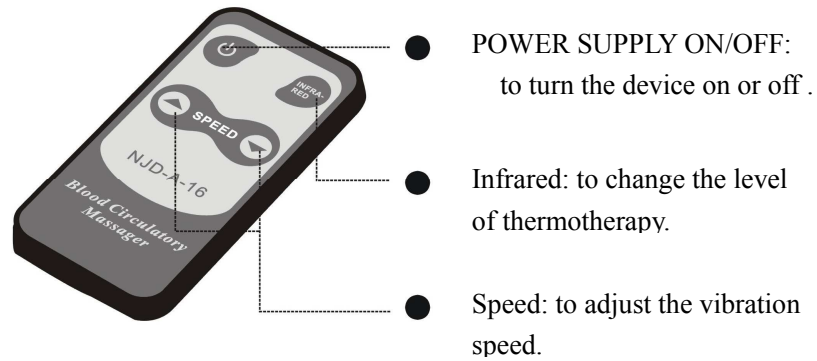
### Use and effectiveness

This device can be used for massaging the feet and other parts of the body. When used for massaging the calves, thigh, back and waist, the vibrating massage can alleviate the feeling of tiredness and provide relief for stiff muscles.

#### ●Foot massage

Sit down in a comfortable sitting position and then rest your feet on the plate. You do not need to remove your trousers or your socks. Feet can be an indication of the

general state of the body's health. Foot massages can thus reduce deposits clearing them away.

## **Use and effectiveness**

### **•Calf massage**

Sit down or lie down in a comfortable position, then place your calf in the centre of the device.

### **•Thigh and buttocks' massage**

Sit down in a comfortable position, and then place your thigh in the centre of the device. Massaging these parts relaxes tense muscles.

### **•Back and waist massage**

Place the massage system on the sofa in a vertical position using the back to support its base. Sit down and apply the section of your back and waist that you wish to massage. If you experience any discomfort, place a towel or thick piece of cloth between the massage device and this section of your back and waist.

### **•Palm massage**

Stimulating certain points of the palm can relieve stiffness felt in the back of the nape of the neck and the shoulders.

#### SUGGESTION FOR USE

In order to relieve certain pain, we recommend following the programmes below:

Health problem	Vibration level		Infrared level		Magnetic wave
	Speed	Time (min)	Intensity	Time	
Joint pain	5	12	3	6	Yes
Sleeping problems	6	15	3	8	Yes
Indigestion	7	15	1	7	Yes
Chills	5	15	1	7	Yes
Head aches	5	15	1	5	Yes

#### **MAINTENANCE AND STORAGE**

Clean the device before any prolonged period of time during which it will not be used. Please turn the power supply switch to the off position after each use of the

TRI MAX.

## **CLEANING**

- Ensure that the switch is in the stop position and that the plug has been pulled out before cleaning the device.

## **MAINTENANCE AND STORAGE**

- Do not use formaldehyde or organic solvents for cleaning the device, so as to avoid any change in the leather and any cracks in the plastic.
- In case of stubborn oil stains, clean using a cotton cloth with detergent followed by a dry cloth for wiping. The cleaning should be done using a clean cloth.
- A damp cloth cannot be used for cleaning the remote control and the area around the switch. Please use a soft dry cloth.
- Avoid any detergent entering into the ventilation window located in the back.
- Do not wash the unit with a cleaning brush or hard object as you may damage the unit.
- Refrain from disassembling the device to avoid any risk.

## STORAGE

Unplug the device from the plug socket and let it cool down before storing it in its case or in a clean dry place, beyond the reach of children. Do not roll the power supply cord around the unit, as the cable may damage and break.

## TROUBLE SHOOTING

Malfunction	Cause	Solution
Does not work after being switched on	1. The power supply does not work. 2. The power supply cable may be damaged.	1. Check if there is any electric current. 2. Send the device back to the technical centre
	3. The button does not work	3. Switch to ON/OFF: the device will work after a "twitch" is heard.
Stops after a	The operating time has	Turn the unit on again

while	expired.	
-------	----------	--

### **Technical data**

- Product: BLOOD CIRCULATORY MASSAGER
- Model : NJD-A-16
- Input : 220V~ 50Hz
- Power : 60W
- N.W. : 3.8Kg

### RECOMMENDATIONS CONCERNING ELIMINATION OF WASTE

- The packaging consists exclusively of materials that are harmless for the environment and can be left with your local waste sorting centre to be used as recycled materials. The cardboard can be put in the paper collection container. The packaging films should be given to local waste sorting centre.

- When you no longer intend to use the appliance,  
discard it in an environment-friendly way, in keeping  
with local regulations.

**SEVEN SPORT, s. r. o.**

Bořivojova 35/878 130 00 Praha 3, ČR IČO: 268 47 264, DIČ: CZ26847264  
Objednávky: +420 556 300 970, objednavky@insportline.cz  
Reklamacie: +420 556 770 190, mobil: +420 604 853 019, reklamacie@insportline.cz  
Servis: +420 556 770 190, mobil: +420 604 853 019, servis@insportline.cz  
Fax: +420 556 770 192, (servis +420 556 770 191)  
Web: www.insportline.cz, www.worker.cz, www.worker-moto.cz

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**INSPORTLINE, s.r.o.**, Bratislavská 36, 911 05 Trenčín, IČO: 36311723, DIČ: SK2020177082  
Objednávky: +421(0)326 526 701, +421(0)917 649 192, objednavky@insportline.sk  
Reklamacie: +421(0)326 526 701, +421(0)918 408 519, reklamacie@insportline.sk  
Fax: +421(0)326 526 705  
Web: www.insportline.sk, www.worker.sk, www.worker-moto.sk